

Advanced AED and multi-functions

Defi 5S



Lead the way
to save a **life**



Advantages

- 2 Languages to choose**
- 5 Levels of sound volume**
- Adult / Child Mode Button**
- Infra-Red for data transmission**
- Daily, Weekly and Monthly Self-test**
- Using the AED while in the bag**

How Easy to use

Start the device

Turn on the power.

The voice instruction starts immediately for guide on how to operate AED.



STEP 1

Attach the pads to the patient

After attaching the pads, AED will immediately start analysing ECG.

The device will move to next step only after you attach the pads well.



STEP 2

Select energy

International resuscitation guideline* recommends to reduce the energy for child patient to avoid burning heart muscle.

Meditech Defi 5S series can reduce the energy to 30/70J just by select the energy range.



STEP 3

Press the button

If the shock is needed, the shock button is flashing. Press the button to deliver the electrical shock to the patient.



STEP 4

[illegible]

Print page can preview selected page in the playback interface, which can be used for printing basis information.




Water resistant buttons and soft



Technical Specifications

Defi 5S

Self-test Mode	Daily, Weekly, Monthly	Power Type	Battery non-rechargeable
Waveform	Adult, Child	Dimensions	12V, 2800mAh Li-MnO ₂ Cell
Energy	Biphasic truncated exponential	Weight (with battery)	Depth (D): 80mm
Energy Sequence	200-Joules max.	Operating Temperature	Width (W): 240mm
	Child Mode:	Operating Humidity	Length (L): 300mm
	50/50/75 Joules		1.9Kg (4.2 Pounds)
	Adult Mode:		0°C to 40°C
	150/150/200 Joules		Relative humidity between 30% and 95% (non-condensing)
Charge time (New, at 25°C)	Less than 6 sec. to 150J		-20°C to 55°C
Voice Prompt	Less than 8 sec. to 200J		
Visual Indicators Control	Extensive voice prompt		
ECG storage	LED prompts		
Data transmission	Two buttons: ON/OFF, Shock		
	1,500 events.		
	Infrared		

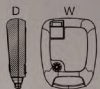

Accessories

Carry Case



Battery

CR123A - 4x2 2.8Ah 12Volt non-rechargeable Battery



Disposable Pad

Designed for the use of external defibrillation, non-invasive pacing, synchronized cardioversion and monitoring on adults and children weighing more than 20 kg (55 lbs).

